

Perfecting Performance



Occlutech® **Delivery Set**




Product Range - Dimensions and sizes.

Art No.	Size	Length Sheath	Length Loader	Dilator Tip	Hemo-valve	Distal Angle
51DS007	7F	800mm	120mm	25mm	70mm	45°
51DS008	8F	800mm	160mm	25mm	70mm	45°
51DS009	9F	800mm	160mm	25mm	70mm	45°
51DS010	10F	800mm	160mm	25mm	70mm	45°
51DS011	11F	800mm	200mm	25mm	70mm	45°
51DS012	12F	800mm	200mm	25mm	70mm	45°

Perfecting Performance

- The micro-braided internal support structure gives optimal stability and excellent torquability.
- Inner surface designed to optimize performance.
- Kink resistant.
- Radio-opaque.

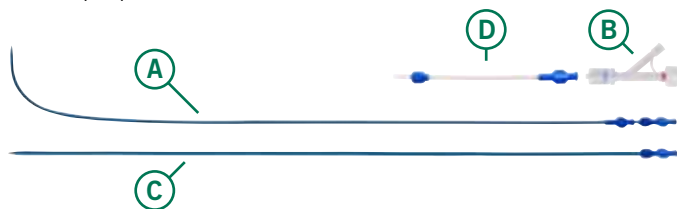

Occlutech delivery set components.

- (A) Delivery sheath
- (B) Hemostasis valve
- (C) Dilator
- (D) Transparent loader

Sheath

The sheath (A) is of Mullins type and consists of a smooth outer layer and a low friction inner layer reinforced with metal braiding. The braid improves the strength of the sheath and provides radiopacity.

The sheath has a proximal female connector, for connection to the dilator (C) or to the loader (D).



Dilator

The dilator (C) is made from smooth conformable low friction radio-opaque material. The dilator is inserted into the sheath (A), to form a supporting inner body.

The inner lumen of the dilator has 0,035" guide wire compatibility.

Loader

The loader (D) is made of low friction material and transparent to allow inspection for air bubbles after pulling the occluder into the loader.



Transparent loader.

Customer Service International

Occlutech, Rönnowsgatan 8, S-252 25 Helsingborg, Sweden
Tel +46 42 33 65 25, Fax +46 42 311 09 70
order@occlutech.com, www.occlutech.com